

Technical Specifications

General Information

Throughput Up to 180 tests/hour

Measurement Principle

Micron superparamagnetic particles platform
With alkaline phosphatase (ALP) labeled reagents
And AMPPD substrate

Sample Assembly

Sample Feeder Continuous sample loading and offloading

Sample positions 50

Sample Probe Steel probe, liquid level detection, clot detection,
Horizontal and vertical collision protection

Sample Volume 10~200ul, with increment of 1ul

Sample Dilution Automatic dilution with 1:2~1:40

Measurement and Reaction System

Detection Mode	Photon counting
Signal Detector	Photomultiplier (PMT) photo counter
PMT Calibration	LED reference module
Incubation positions	82
Temperature	Accuracy: 37±0.3 °C
Mixing Unit	Non-touch vortexer mixing

Reagent Assembly

Reagent Disk 15 reagent positions, continuous loading,
real-time mixing for magnetic particles

Cooling System 2~8 °C

Reagent Pack 50 tests/pack, 100 tests/pack

Reagent Scanning Built-in/External bar code reader

Reagent Probe Steel probe, liquid level detection, clot detection,
Horizontal and vertical collision protection

Reagent Volume 10~200ul, with increment of 1ul

Working conditions

Power Supply	100-240V~50/60Hz
Input Power	500VA
Operating Temperature	15~30 °C
Relative Humidity	35%RH~85%RH, without condensation
Dimension	860mm*740mm*560mm (W*D*H)
Weight	132Kg
Altitude Height	-400m~3000m (70kpa)

Substrate System

Substrate Bottle 300/500 tests per bottle, bottom piercing

Substrate Loading Capable of continuous loading during testing

Substrate Heating Constant heating

Aspiration Volume 200ul

Magnetic separation unit

Separation Assembly 3-phase magnetic separation

Temperature Accuracy: 37±0.3 °C



CL-960i

Chemiluminescence analyzer

180 Test/hr



'Precise, Fast & Compact'
Fully automated CLIA analyzer

www.rapiddiagnostic.in

For Distributor, Wholesaler & Retailer Enquiries, contact:

Aspen Diagnostics Pvt. Ltd.

(A Rapid Diagnostic Group of Company)

B - 82, G. T. Karnal Road Industrial Area,

Delhi - 110 033, INDIA.

T: +91-11-45201100 | F: +91-11-27413437


drjindal@rdgc.in | www.rapiddiagnostic.in

Local Contact:


CL-960i The right size for your lab



Precise magnetic separation with consistent performance reliability

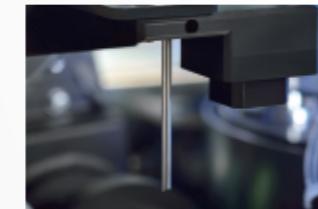


Solid-liquid waste separation
Environmental friendly

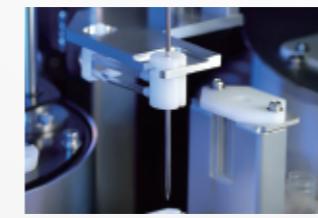

Wiper washing technology
Enhance efficiency and minimize carryover



Continuously loading of
reagents and consumables

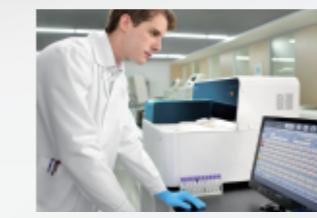

User maintenance free
Lighten the workload



One-key start
Extremely easy to use



Intuitive software interface
Easy access to all functions



Intelligent consumables management

Test Menu

Thyroid	Infectious Disease	Fertility	Tumor marker	ToRCH	Hypertension	Anemia
FT3	HIV Combo*	Total β-HCG	CEA	Toxo IgG*	Angiotension I*	Ferritin
FT4	HBsAg*	FSH	AFP	Toxo IgM*	Angiotension II*	Vitamin B12
T3	Anti-HBs*	LH	CA125	Rubella IgG*	Renin*	Folate
T4	HBeAg*	Prolactin	CA15-3	Rubella IgM*	Aldosterone*	RBC Folate
TSH	Anti-HBe*	E2	CA19-9	CMV IgG*	DS' Screening	Growth Hormone
Anti-TG	Anti-HBc*	E3	FPSA	CMV IgM*	PAPP-A*	hGH*
Anti-TPO	Anti-HCV*	TESTO	TPSA	HSV-1/2 IgG*	Free β-HCG*	IGF-1*
Tg	Anti-TP (Syphilis)*	PROG	NSE	HSV-1/2 IgM*		
rT3	HBeAg quantitative*	Free testosterone*	CYFRA 21-1	HSV-1 IgG*	Liver Fibrosis	
TRAB*	Anti-HBc IgM*	17-OH PROG*	CA72-4	HSV-2 IgG*	LN	
		SHBG*	PG I*		HA	
Bone	Cardiac	AMH*	PG II*		PIIINP	
Metabolic	CK-MB	SCCA*	HE4*		CIV	
Intact PTH	Myoglobin	ProGRP*	CA50			
Calcitonin	Troponin I	Insulin	CA242		Inflammation	
Vitamin D total	BNP	C-peptide			PCT	

* Coming Soon